


iFIX ♦ GENERAL PRODUCT INFORMATION

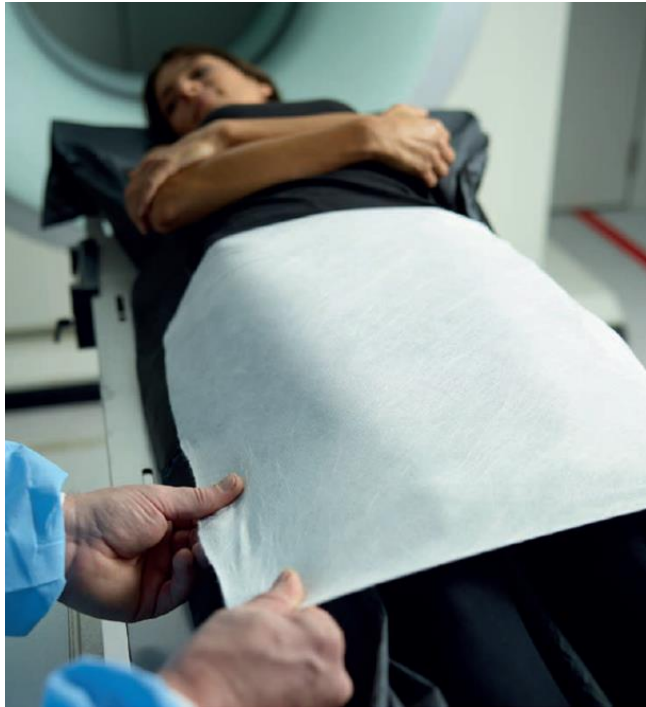

iFIX is a patient stabilization system comprised of both reusable and disposable components.

The comfortable stretchy fleece is intended for single patient use.

It contours to the patient’s anatomy like a second skin and secures them in the desired position.

The iFIX Fleece as a single use component improves hygienic conditions and lowers the risk of cross-contamination.

The iFIX Fleece is attached to a reusable anchor product.


Depending on your table’s design, you can choose between the iFIX Patch that adheres to the table surface, the iFIX Slot Adapter that slides into the channels, or the iFIX Patch Adapter which connect to the fleece products.

The unique micro anchors of the iFIX Patches and Adapters strongly secure the iFIX Fleece to ensure optimal patient stabilization thus minimizing motion artifact.

iFIX PRODUCT LINE

disposable fleece products

reusable anchor products

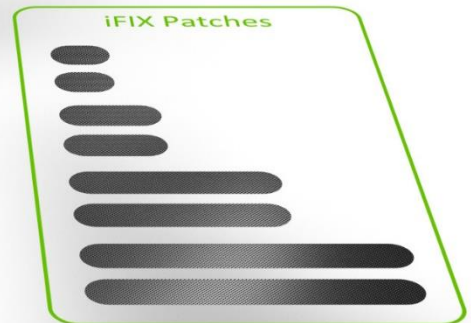

iFIX ♦ PRODUCT PORTFOLIO

FLEECE PRODUCTS (SINGLE PATIENT USE)



iFIX Fleece 15 consists of skin-friendly stretchy fleece material and comes from a role with 15 cm width. The textured fleece side adheres to the iFIX anchor products. The fleece contours to the patient's anatomy like a second skin and stabilizes them in the desired position.



iFIX Fleece 50 consists of skin-friendly stretchy fleece material and comes from a role with 50 cm width. The textured fleece side adheres to the iFIX anchor products. The fleece contours to the patient's anatomy like a second skin and stabilizes them in the desired position.



iFIX Head Mask consists of skin-friendly stretchy fleece material and come in different sizes (**XS, S, M, L, XL**).

The textured fleece side adheres to the iFIX anchor products. The fleece contours to the patient's head like a second skin and stabilizes the head in the desired position.

ANCHOR PRODUCTS (REUSABLE)



iFIX Patch serves as anchor for the iFIX Fleece products. The side with the adhesive tape adheres to the table surface. Unique micro anchors on the black surface strongly secure the iFIX Fleece to ensure optimal patient stabilization.

iFIX Patches are available in **15, 25, 50** and **75** cm.



iFIX Slot Adapter serves as anchor for the iFIX fleece products. The iFIX Slot Adapter slides into the channels of tables. Unique micro anchors on the black surface strongly secure the fleece to ensure optimal patient stabilization.

iFIX Slot Adapters for **4.0** and **9.5** mm channels are available.



iFIX Patch Adapter serves as anchor for the iFIX Fleece products. The black loop side of the iFIX Patch Adapter adheres to unique micro anchors or to the iFIX Patches.

Unique micro anchors on the grey surface strongly secure the fleece to ensure optimal patient stabilization.